JAY[®]

EASY VISCO CUSHION

JAY[®] Easy Visco Cushion

Important Consumer Information

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LT 73

SI 77

NOTE: Check all parts from shipping damage. In case of damage DO NOT use. Contact Carrier/Sunrise for further instructions.

Sunrise Medical recommends that a clinician such as a doctor or therapist experienced in seating and positioning be consulted to determine if a JAY Easy Visco Cushion is appropriate. Cushions should only be installed by an authorized Sunrise Medical supplier.

⚠ GENERAL WARNINGS!

- DO NOT install this equipment without first reading and understanding this manual.
- Keep this manual with the product for future reference.
- If you are unable to understand the Warnings, Cautions and Instructions, contact a qualified clinician or supplier - otherwise injury or damage may occur.
- Hook and Loop velcro should not come into prolonged contact with moisture as this may degrade the adhesive and lead to a failure of the velcro.
- Installing a cushion on a wheelchair may affect the center of gravity of the wheelchair and may cause the wheelchair to tip backwards, potentially resulting in injury. Always assess for the potential need for antitippers or amputee axle adapter brackets to be added to the wheelchair to help increase stability.
- Prior to prolonged sitting, any cushion should be tried for a few hours at a time while a clinician inspects your skin to ensure that red pressure spots are not developing. You should regularly check for skin redness. The clinical indicator for tissue breakdown is skin redness. If your skin develops redness, discontinue the use of the cushion immediately and see your doctor or therapist.
- The JAY Easy Visco Cushion is designed sitting comfort and postural support, remember no cushion can completely eliminate sitting pressure or prevent pressure sores.
- If you need more pressure distribution for skin protection, ask your authorised supplier for a free demonstration of any of the other Jay cushion products.
- The JAY Easy Visco Cushion is not a substitute for good skin care including; proper diet, cleanliness, and regular pressure relief.
- Allow the cushion to warm to room temperature before using. Sitting on extremely cold or hot surfaces can cause skin damage.
- Avoid sharp objects or exposure to excessive heat or open flame.
- OBSTRUCTIONS: Do not place any obstructions between the user and the cushion because it will reduce product effectiveness.
- CUSHION AND COVER ORIENTATION: Product must be used with the cushion cover facing up. If the cover is not used correctly it may reduce or eliminate the cushion's benefits and could increase risk to the skin and soft tissue.

JAY EASY VISCO

The Easy Visco Cushion is a wheelchair cushion that aids in the stabilization of the pelvis and torso of the wheelchair user and properly positions the wheelchair user and properly positions the wheelchair user's body while seated.

The Easy Visco cushion is designed for high level comfort and positioning with mild pressure redistribution. It incorporates a lightweight pre-contoured foam with media/lateral thigh support and Visco elastic foam in the seat well to allow immersion of the ischia's (seat bones) and load distribution to the femurs. The area of the sacral is notched. The outer cover is moisture resistant with a no-slip bottom to keeping the cushion in place during daily use and transfers.

The intended lifetime of this product is 5 years.

Maximum user weight for this product is:

150 kg (330 lbs)

FITTING THE JAY EASY VISCO CUSHION

- Select a Jay Easy Visco Cushion that matches your hip width.
- The cushion length should extend to within 25 mm - 50 mm, (1" - 2"), from the back of the knees.
- Place the cushion to the back of the chair with the material identification tags and Visco elastic foam towards the rear.
- When fitted correctly, your hips should be to the back of the chair.
- Your ischials (seat bones) should be centered on the visco elastic foam in the seat well of the cushion.
- Proper footrest adjustment may enhance sitting comfort and help lower peak sitting pressures.
- Adjust the footrests so your legs rest firmly but comfortably on the cushion.

NOTE:

- Do not use the wheelchair's seat width as a guide in determining the appropriate cushion width. Improper cushion fit may be detrimental to the performance of a pressure distributing cushion and its support of the body.
- Avoid using the cushion on overly stretched seat upholstery. Stretched upholstery does not provide adequate support for the cushion and improper positioning may result. Replace stretched upholstery or purchase the Jay adjustable solid seat.
- Avoid using overly stretched backrest upholstery. Stretched back upholstery does not provide adequate support for the body and improper positioning may result. Replace stretched backrest upholstery or purchase a back support product from Sunrise Medical that is appropriate for you. Contact an authorized supplier for a free demonstration of these products.

EASY MAINTENANCE & CLEANING

Maintenance

Monthly cleaning and regular maintenance may help extend the life of your cushion. During cleaning, component inspection is recommended. Check the cover for tears and excessive wear. Check the foam base and the Visco foam in the seat-well to ensure foam consistency.

Washing

TO CLEAN COVER FIG.1.0

- Unzip the cover, “A” and Remove the cover from the base.
- Turn the cover inside out and machine wash in warm water and drip or tumble dry on low heat (follow instructions from cleaning chart using common laundry detergent).
- No ironing required.

Cleaning

- Remove the cover from the foam base
- Spray the surface with a 70% Alcohol based or 0.1% Chlorine solution
- Let it dry
- Wipe off the solution completely with a moistured non abrasive cloth with clear water
- Let it completely airdry

⚠ CAUTION!

- Do not dry clean covers or use industrial washers and dryers to clean cover.
- Do not steam autoclave.
- Do not bleach.

After washing, fit cover back over foam base. Ensure the back of the cover matches up to the back of the foam base.

NOTE:

- If utilized in an institutional setting, write the resident's name with a permanent marker on the content label of the covers. This will help to ensure the return of the covers after washing.
- Use only Jay-designed covers in conjunction with Jay Easy Visco Cushions.

TO CLEAN FOAM BASE FIG.1.0 - 2.0

- Remove the cover from the foam base.
- Wipe with warm water and soap.
- Wipe off with a clean cloth.
- Avoid complete submersion of the foam base in water.

NOTE: See Detailed Cleaning Instructions overleaf.

HYGIENE MEASURES WHEN BEING RE-USED

This JAY cushion is designed to be re-used. The frequency of re-use depends largely on the way the cushion has been used previously. It is essential that the cushion is properly cleaned before re-use according to the local authorities regulation (see section “Easy Maintenance & Cleaning”). It is strongly recommended to replace the cover when reusing the cushion.

CONDITIONS FOR SAFE STORAGE

- Products should be stored in a dry place: at 20 – 75% Relative Humidity.
- Protect from direct sunlight and dust for long periods.
- Away from direct heat and at a temperature between 5°C (40°F) to 40°C (104°F)
- Position and stack the packaged products according to the “this side up” arrows on the packaging.
- If a package is high and narrow and therefore potentially unstable, ensure it is well secured on the pallet or in the racking, to prevent it from falling down.
- Do not stack anything on top of the packaged products when the “do not stack” mark on the packaging is shown.
- Product should never be stored outside and subject to the elements.

JAY EASY VISCO CUSHION WARRANTY

Each Jay cushion is carefully inspected and tested to provide peak performance. Every Jay cushion is guaranteed to be free from defects in materials and workmanship for a period of 24 months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within 24 months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears or burns, nor to the cushion's removable cover.

The removable cover is guaranteed to be free from defects in materials and workmanship for a period of six (6) months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within six (6) months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge.

Claims and repairs should be processed through the nearest authorized supplier. Except for express warranties made herein, all other warranties, including implied warranties of merchantability and warranties of fitness for a particular purpose are excluded.

Model	JAY Easy Visco
Article Number*	JEVwwddF (Flat), JEVwwddC (Contour)
Weight in ready to use state (Kg)**	2.7
Anti Slip Coating	Yes, at the bottom of the Outer Cover
Outer Cover Material	Polyester
Insert Material	N/A
Inner Cover Material	N/A
Foam Base Material	Polyurethane (PU) Foam
Foam Density	74.5 Kg/m3
Compression Hardness	21 +/-3lbs (9.5 +/-1.36 Kg), 50% Deflection with 50 mm
User Min Weight (Kg)	No minimum weight
User Max Weight (Kg)	150
Spatial Dimension in ready to use state*	Seat Widths 250 to 600 mm, Seat Depths 250 to 600 mm, Height approx 60 mm

*For the complete offered sizes (Seat Widths - Seat Depths) and different components, please refer to the Order Form of the product. The scope of product delivery will vary according to the clinical assessment and prescription made from a clinician.

**Data based on an average configuration

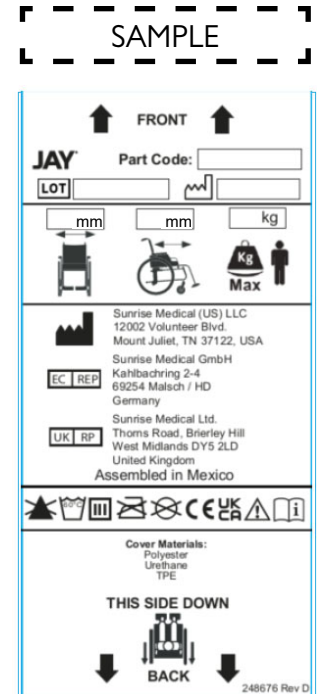
N/A: (Not Applicable) as the cushion model does not include the element.

ISO 7010-M002 -
Instruction manual/
booklet must be read!
(Blue Icon)

Type:	Product Name/ SKU Number	XXXX-XX-XX	Date of Manufacture
LOT	LOT: defined by Julian Date	Part Code	Assembly Part Number
XXX mm	Seat width	MD	This symbol means Medical Device
XXX mm	Seat depth		Manufacturer's address
XXX kg	Max user weight	EC REP	European Authorised Representative
CE	CE Mark	UK RP	UK Responsible Person
UK CA	UKCA Mark	CH REP	Swiss Representative's address
	Consult instructions for use		Importer's address

The management system of SUNRISE MEDICAL is certified to EN ISO 13485 and ISO 14001.

We at SUNRISE MEDICAL have been awarded the ISO-13485 certificate, which affirms the quality of our products at every stage, from R & D to production. This product meets the requirements in accordance with EU and UK regulations. Options or accessories shown are available at extra cost.

Indications

The varieties of cushions fitting variants, as well as the different positioning components to deliver higher positioning and skin protection benefits according to each user requirement and indication, mean that it can be used by a range of users with different positioning and skin protection needs and indications as follows [R7-R14] :

- Spinal Cord Injury
- MS
- Neurological disease pattern (hemiplegia, M.Parkinson)
- Extremity defect deformity
- Joint contracture
- Amputation
- Tetraplegia

Contraindications

The cushion for a wheelchair shall not be used in case of [R3, R4]:

- Redness developed in the skin contact area
- Decubitus

NOTE: General user advice. Not following these instructions may result in physical injury, damage to the product or damage to the environment!

Notice to the user and/or patient: any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

As the manufacturer, SUNRISE MEDICAL, declares that this product conforms to the Medical Device Regulation (2017/745).

As the manufacturer, SUNRISE MEDICAL, declares that the product conforms to the UK Medical Devices Regulation 2002 No. 618.

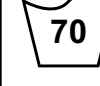
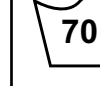
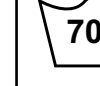
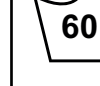
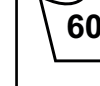
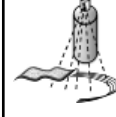
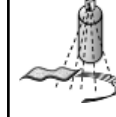
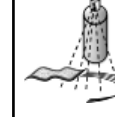
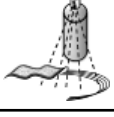
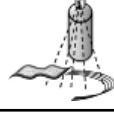
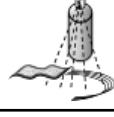
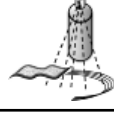
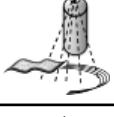
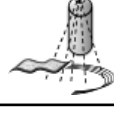
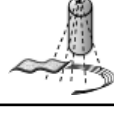
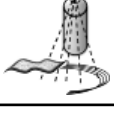
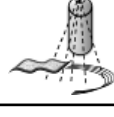
B4Me special adaptations

Sunrise Medical strongly recommends that in order to ensure that your B4Me product operates, and performs as intended by the manufacturer; all the user information supplied with your B4Me product is read and understood, before the product is first used. Sunrise Medical also recommends that the user information is not discarded after reading it, but it is kept safely stored for future reference.

Medical Device Combinations

It may be possible to combine this Medical device with one or more other Medical Device or other product. Information on which combinations are possible can be found at www.SunriseMedical.co.uk. All combinations listed have been validated to meet the General Safety and Performance Requirements, section 14.1 of the Medical Device Regulation 2017/745.

Guidance on the combination, such as mounting, can be found at www.SunriseMedical.co.uk.

Approved Cleaning Procedures for JAY Cushions									
Product	JAY Lite	JAY J3	J2 - J2 DC - GS	JAY Easy Visco	JAY Soft Combi P	JAY Xtreme Active	JAY Basic	JAY Easy	JAY Balance
Outer Cover (Microclimatic and Stretch)									
(1). Wet Wash 70°C 									
Outer Cover (Incontinent)									
(2). Wet Wash 60°C 									
Inner Cover									
(3). Spray and Wipe disinfection  or Wet Wash 60°C	NA	NA	NA	NA	NA		NA	NA	
Cushion core (foam base)									
(4). Steam autoclave up to 105 °C (1) 	✗	✗	✗	✗	✗				✗
Spray and wipe disinfection 						✗	✗	✗	
(5). Immersion disinfection 				✗	✗	✗	✗	✗	✗
Cushion elements									
Air element	NA		NA	NA	NA	NA	NA	NA	
Positioning elements	NA			NA	NA		NA	NA	
Fluid pad	NA			NA	NA		NA		
1. 70°C is the maximum recommended temperature. 2. 60°C is the maximum recommended temperature. 3. The Xtreme Active and Balance inner covers can be washed at 60°C as well as spray and wipe disinfected. Possible discoloration if washed at 60°C but the functional aspect of the covers is not diminished. 4. An autoclave is a self locking machine that sterilizes with steam under pressure. 5. At low temperature (< 30°C) 6. NA (Not Applicable) as the cushion model does not include the cushion element									
The product warranty shall not apply to products damaged due to the failure to follow these cleaning instructions or the use of inappropriate cleaning compounds.									
Recommended Disinfectants									
Solutions containing: Alcohol 70% solution base, Chlorine 0,1% solution (must be rinsed always with clear water). Refer to the cleaning instructions procedures included in this manual									

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Sunrise Medical B.V.
Sunrise Medical HCM B.V.
Sunrise Medical S.A.S