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Head Rests and Holders

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Poggiatesta e staffe

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CH-REP. CH-Importeur: INNOVAhelp GmbH, Thunstr. 5a, CH-3110 Münsingen

01. Purpose and Indication

The Rolko head rests and the corresponding holders support the positioning of the head of persons in need in wheelchairs. For this purpose, the holders can be attached to the push handles, the back plate or the seat shell of wheelchairs. This instruction manual is aimed at trained specialists in the field of rehabilitation and care. The necessity and suitability of the head rests and the corresponding holders must always be evaluated by qualified personnel.

02. Contraindication

The head rest holders HA-KO-ALU-C22, HA-KO-ALU-C25, HA-KO-ALU-KG-C22 and KO-ALU-KG-C25 are only suitable for tube diameter of 22 or 25 mm, as the head rest holders HA-UV-KO-S and HA-UV-KO-38-S are suitable for tube diameters from 22 up to 25 mm. The head rest holders are not suitable for all other tube diameters. There must be enough space on the push handle tubes of the wheelchair for the attachment of the head rest holders. The usage of other operating elements, such as brake levers, must not be impaired by the head rest holders. For the head rest holders HA-KO-ALU-R, HA-KO-ALU-KG and U-BV-KOU-S, the back plate or seat shell must be suitable for the attachment. The suitability is to be evaluated by trained specialists. The required fixation material for this is not part of the scope of delivery. These head rests are unsuitable for folding wheelchairs and wheelchairs with a textile backrest. The Rolko headrest holders are not approved for use in motor vehicles.




03. Safety Instructions

- These instructions must be read and followed before installing the product.
- Check the head rest holder and the head rest for damages before using them for the first time.
- Damaged / defective head rest holders or head rests must never be used and need to be replaced immediately.
- The use in a motor vehicle is not permitted.
- The instructions for use are part of the head rest holder and must be passed on.
- When installing a head rest holder on the push handles, the use of other operating elements, such as brake handles, must not be impaired under any circumstances.
- If serious incidents occur in connection with this product, these must be reported to Rolko Kohlgrüber GmbH and to the competent authority of the EU Member State in where the product is being used.

04. Assembly and Adjustment

Head rest holder HA-KO-ALU-C22 and HA-KO-ALU-C25:



The two parts C of the head rest holder are inserted into the tube of part A on both sides and clamped in the correct length with the clamping lever of part A. The fixation clamps of part C are attached to the wheelchair push handles as shown in the illustration below on the right. The part B of the head rest holder is inserted into the middle clamping piece of part A and clamped at the desired height with the help of the clamping lever. The head rest is attached to the ball head with the clamp U-BST-KO-5/6-S and three M5 screws. The head rest can be adjusted in the height via the middle clamping piece of part A and in depth via the two joints on part B. To do this, loosen the two M8 Allen screws and retighten them in the desired position.

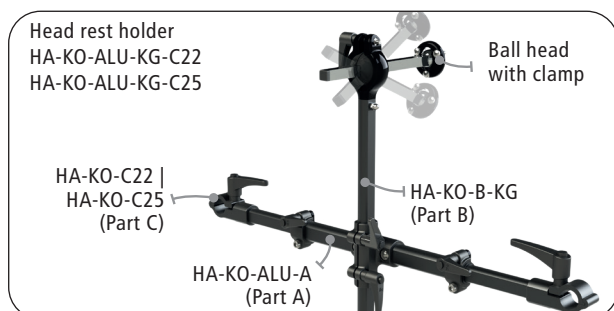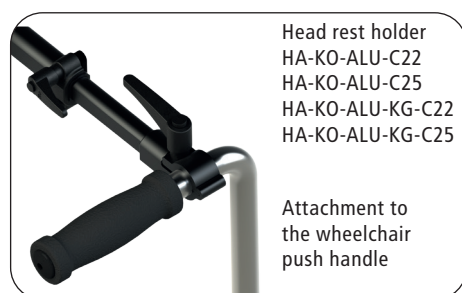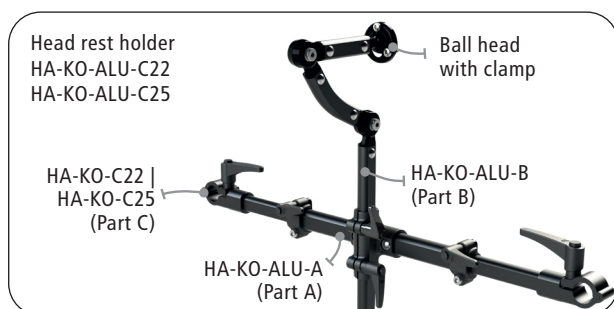Head rest holder HA-KO-ALU-KG-C22 and HA-KO-ALU-KG-C25:
The attachment and adjustment of those head rest holders work as with the models HA-KO-ALU-C22 and HA-KO-ALU-C25. The depth and inclination of the position of the head rest is set here on the ball joint of part B. To do this, loosen the clamping lever on the upper ball joint of part B, position the crossbar at the correct depth and incline and retighten the clamping lever at the end. The head rest is attached to the ball head with the U-BST-KO-5/6-S clamp and three M5 screws.




Head rest holder HA-UV-KO-S and HA-UV-KO-38-S:

The attachment to the push handles and the adjustment work similarly to the head rest holders HA-KO-ALU-C22 and HA-KO-ALU-C25. Instead of the clamps for attachment to the push handles, there are clips here, which are tightened to the push handles with a star grip. The clips are suitable for tube diameters from 22 to 25 mm. The height and depth of the head rest can be adjusted by loosening the star handles on the upper vertical tube, moving the vertical tube and the short upper tube with the ball mount into the desired position and fixing them by tightening the star handles. With HA-UV-KO-38-S, the attachment to the push handles is solved differently: Here the tube openings from part C are placed over the push handles and then fixed behind the grip surfaces of the push handles by the star handles.

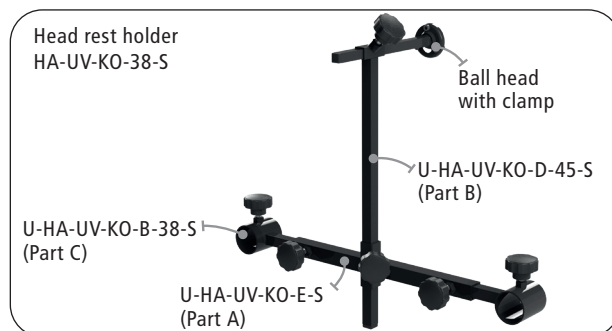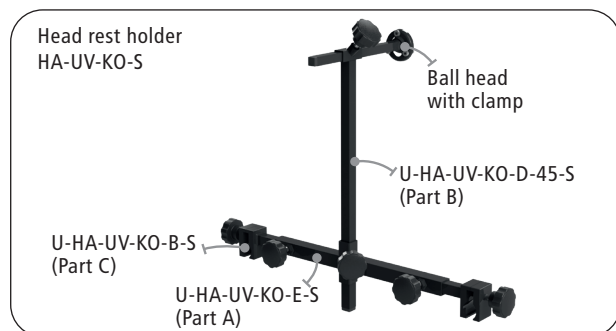Head rest holder HA-KO-ALU-R, HA-KO-ALU-KG and U-BV-KOU-S:

The mounting plate with the clamp for the vertical tube of part B is firmly screwed to the back plate or seat shell, depending on the type of the back plate or seat shell. The fixation material is not included in the scope of delivery. The attachment and adjustment of the height and depth of the head rest are carried out in the same way as with HA-KO-ALU-C22 and HA-KO-ALU-C25.

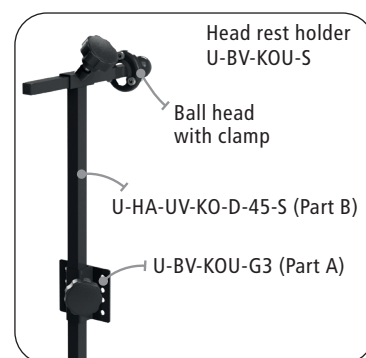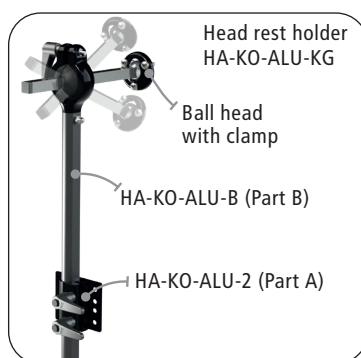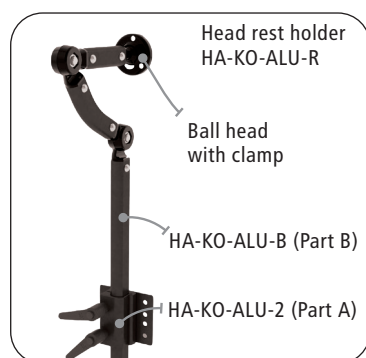Head rest models

05. Cleaning and Disinfecting

For cleaning, the head rest holder can be wiped off with a damp cloth with the addition of a mild detergent, if required. Use lukewarm water and, if necessary, non-abrasive, pH-value neutral, commercially available household cleaners. If necessary, the head rest holder can be cleaned with common household disinfectants by spray or wipe disinfection. The following cleaning agents and disinfectants or processes that contain the following are not permitted for the polyurethane head rests:

- hot water, steam
- concentrated acids and alkalis
- alcohol
- chlorinated hydrocarbon
- ammonia

06. Maintenance

Verify the function of the head rest holders before each use and make sure that they are undamaged and in perfect condition.

Before each use, check:

- the tight fit of all fixations
- the correct adjustment of the head rest holder
- the product visually for any damages

07. Re-use

This product is suitable for re-use. Clean the head rest holder according to the specifications under "05. Cleaning and disinfection".

The following points are to be checked:

- if the function of all fixations is given
- the head rest holder and the head rest are visually undamaged
- if all individual parts and the instructions for use are available
- the head rest holder is to be stored protected from dust and moisture

08. Disposal

All parts of the head rest holder and also the head rest must be disposed of safely with household waste. If possible, dismantle the head rest holder into its individual parts to enable the materials to be sorted.


09. Technical Data

Scope of delivery: head rest holder, instruction manual, optional: head rest

Materials used: aluminum, steel, PA6, PU

Temperature range: -20° C to +50° C

Relative humidity: max. 80%

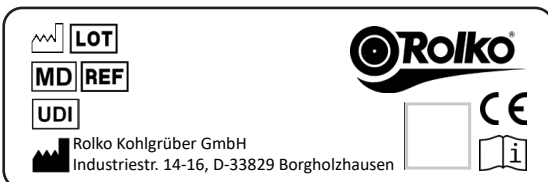 Manufacturer

 Follow the instructions

 Production date

 Item no.

 Batch

 Medical device

 Unique device identification

As an accessory for medical products of class 1, this product has a CE mark in accordance with regulation (EU) 2017/745 (MDR) and thus complies with all German and European safety requirements, that apply to the product. The product fulfils the corresponding sections of standard specification DIN EN ISO 21856. The EU declaration of conformity is available for download at www.rolko.com.

Instruction manual. Head rests: PU-KO-5-S | PU-KO-6-S | PU-KO-7-S | PU-KO-10-S

Head rest holders: HA-KO-ALU-C22, HA-KO-ALU-C25 | HA-KO-ALU-KG-C22, HA-KO-ALU-KG-C25 | HA-UV-KO-S | HA-UV-KO-38-S | HA-KO-ALU-R | HA-KO-ALU-KG | U-BV-KOU-S

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Reprint, also in extracts, only with prior permission by the editor. All rights reserved. Alterations for reasons of technical development reserved. All technical data are characteristic and subject to construction and manufacturing conditioned tolerances. Errors excepted.

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© Rolko Kohlgrüber GmbH, Industriestr. 14 - 16, 33829 Borgholzhausen