Driver

Aluminium Rollator Driver

User's Manual (EN)

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1. General Information

Dear Customer,

Thank you for purchasing our medical device. In this user's manual, you will find all the important information related to the proper and safe use of this product, hereinafter referred to as the "rollator." Before using the rollator for the first time, please read and carefully familiarise yourself with the contents of this manual and follow all our recommendations.

It is especially important to read and follow the safety requirements.

Rehasense continuously improves its products and reserves the right to make changes to the product specifications without prior notice.

2. Intended Use, Indications, and Contraindications

Indications: The walker has been designed as a mobility aid for individuals who experience difficulties with standing and walking. It enables such individuals to return to many everyday activities.

It is intended for people who, due to various medical conditions, require stable support and assistance in moving around. The walker must not be used for purposes other than providing support for standing, walking, and sitting.

Contraindications: The walker should not be used by individuals who have serious balance problems or significant weakness of the skeletal muscles that prevents them from standing or walking independently. The use of the walker should be recommended by a physiotherapist or a doctor.

3. Delivery and General Product Description

The walkers are delivered in cardboard packaging in a folded position. The package contains:

- ✓ folded walker
- ✓ crutch holder set
- ✓ user manual

The DRIVER rollator is equipped with 4 wheels. The lightweight aluminium frame allows it to be folded, stored, and transported.

The product offers easy manoeuvrability, stability, and adjustable handle height.


The frame is coated with powder paint. It is also factory-fitted with a seat and a shopping basket.

4. Requirements and Safety


Before using the walker, read the contents of this manual. The manufacturer is not responsible for any damage or accidents resulting from improper use.

- Maximum Capacity / Use
 - The user's weight must not exceed 150 kg.
 - The standard shopping bag can hold up to 10 kg.
 - The walker is intended for use both indoors and outdoors.
 - The walker must not be used or pushed while seated.
 - Before using the seat, the brakes **MUST** be in the locked position.
- Accident Prevention
 - Use the device only as intended.
 - Use the device only when it is in proper technical condition.
 - Do not make any modifications to the device unless the manufacturer has given written permission.
 - Exercise special caution when adjusting the handle height to avoid trapping fingers or other parts of the body.
 - Use the device only on flat, smooth surfaces.
 - All wheels should always be in contact with the ground. This ensures proper stability and helps prevent accidents.
 - Always lock the brakes when the walking aid is not in use.

IF ANY SERIOUS INCIDENT OCCURS IN RELATION TO THE DEVICE, PLEASE CONTACT YOUR DEALER OR REHASENSE DIRECTLY (contact details at the end of the user manual). The immediate actions should be taken according to the requirements of Regulation of the European Parliament and of the EU Council 2017/745 on medical devices of April 5, 2017.

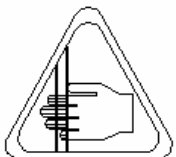 Mind your fingers	Be careful when you see this symbol on the device. Fingers or other body parts of a user and/or assistant can be trapped between the device's parts during folding, unfolding or adjusting of the device. Be extremely careful when adjusting the device. Observe the dimensions of gaps between the device's elements to avoid trapping your fingers or other body parts.
--	---

5. Device Description

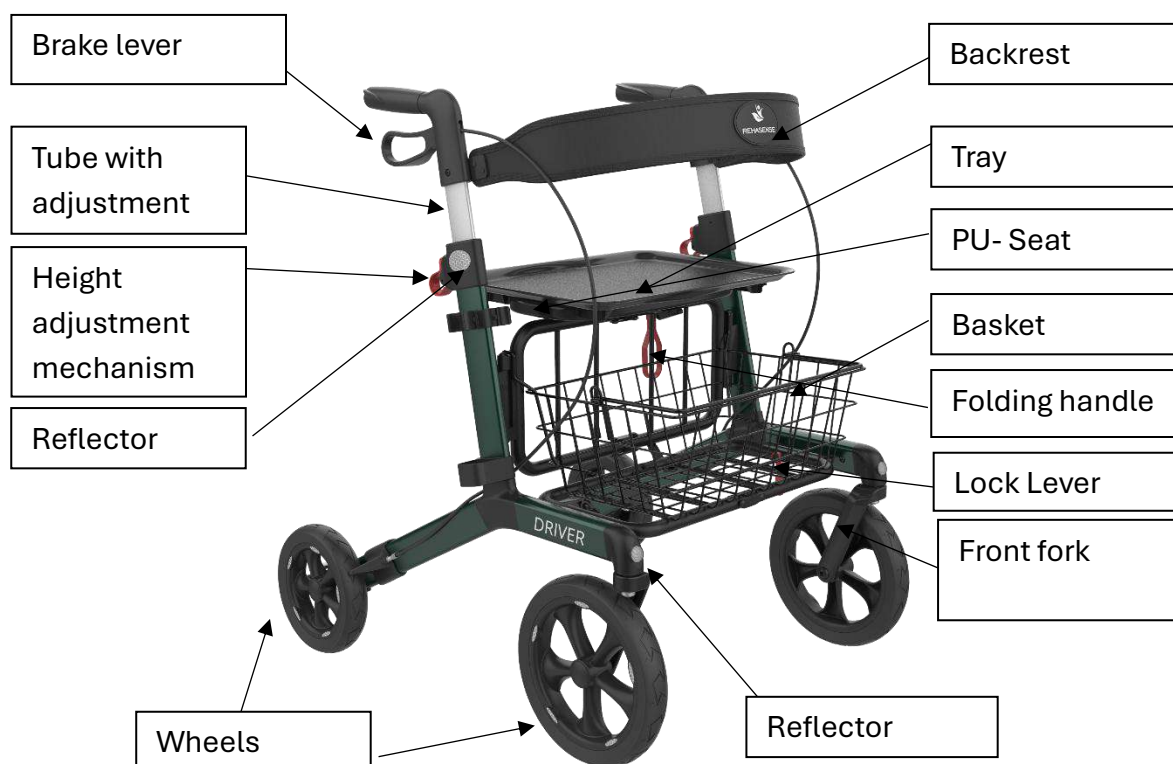6. Use

The above picture shows the DRIVER rollator along with all accessories. After receiving the rollator, follow the instructions below:

Initial Inspection

After removing the device from the packaging, first check that nothing is missing and that all parts shown in the picture are included in the box. If you notice any missing parts, contact the device supplier immediately.

Unfolding the Walker

- Turn the lock lever (visible in the left picture) counterclockwise.
- Pull the handles apart so that the seat rotates clockwise.
- Slide the seat forward until you hear a distinctive “click.”

- The walker is now unfolded and secured against accidental folding.
- See the right picture, which shows the walker in the locked position.
- For safety, check the unfolding by attempting to fold the frame slightly.

- Always ensure, when unfolded, that the strap attached to the seat is secured; it facilitates folding the walker.
- Always check before use that all parts of the walker are in proper working condition.
- Failure to follow these instructions may result in injury to the user.

Folding the Walker

- Pull the strap attached to the seat upward until the frame begins to collapse inward and the red button locks automatically. The walker is now folded.

Adjusting Handle Height

- Adjustment is done by simultaneously pressing and pulling upward the handle height adjustment mechanism.
- Pull one handle, then the other, upward so they are at the same height.
- Both handles should always be set to the same height. Never leave the device with unlocked brakes.
- Correctly positioned handles should be at the user's hip level.

Braking / Releasing the Brakes

- To stop the walker, pull the brake handles upward.
- To lock the wheels for an extended period, press the brake handles downward until you hear a "click."
- Only with the brakes locked may the user turn around and sit on the seat.
- To move again, stand up, turn, and release the brakes by lifting the handles upward.

Safety Rules for Brakes and Seating

- Always keep the brakes locked when the walker is stationary.
- Never sit on the seat without locking the brakes!
- Never lean backward while seated. It is recommended to sit at the front of the seat, taking care not to sit "on" the bag.
- Sit upright and hold the handles.
- Never leave the walker unlocked on slopes or inclines!

Seat Use

- Before use, make sure the brakes are locked.
- Always sit with your back to the direction of movement.

7. Cleaning, Disinfection, and Maintenance

General:

- The walker and its accessories should be cleaned and maintained regularly to ensure optimal performance.
- Light dirt and dust can be removed with a damp cloth.
- Heavier dirt can be cleaned with a mild, warm soap solution.
- Do not use abrasive cleaners or bleach.
- Do not clean with solvents such as gasoline or thinner, as these can damage the paint and plastic components.
- Clean the wheels with warm water and a mild detergent.
- The walker and accessories can be cleaned in a washing cabin at a maximum temperature of 80°C for up to 10 minutes.
- When using high-pressure cleaning, do not direct the jet directly at the wheel bearings.
- The walker and accessories can be disinfected using a soft cloth or sponge soaked in a disinfectant containing 70–80% ethanol. Disinfectants based on chlorine or phenol are not allowed.
- Regularly check that hair or other debris does not accumulate between the wheels and forks.
- Use only disinfectants approved for use (the current list is available on the Robert Koch Institute website: www.rki.de).

Wheels:

- Clean the wheels with warm water and a mild detergent. Dry them with a clean cloth.

8. Sample Label

The product marking in the form of a label is located on the main frame and looks as shown below.

THE LABEL IS FOR REFERENCE ONLY — THE CORRECT ONE IS LOCATED ON THE ROLLATOR FRAME.

	Medical Device	Familiarise yourself with UM	
Label ID	Label id: 5410QQ rev 2	CE MDD	Production Date
Product Name (size)	Driver (L)	2025-08-26	Item number
Barcode EAN 13 in standard GS1	5 907467 824260	REF RRDRVDGL40	Max. User Weight
Barcode standard GS1 (UDI-DI)	UDI (01)05907467824260(11)250820(21)0001	150 kg	Basest max. load
	Can be used outdoor	10 kg	Product Width
		66.8 cm	
		Manufacturer	

9. Technical Data

“DRIVER”	
Maximum user weight	150 kg
Maximum bag load	10 kg
Total weight without bag	9,5 kg
Basket weight	0.948 kg
Length	79 cm
Height (unfolded)	80 cm
Width (folded)	22 cm
Width (unfolded)	67 cm
Seat width	48 cm
Handle height (adjustable)	80 - 97 cm
Seat height	59 cm
Brakes	Brakes with a rear wheel locking function
Wheels	Soft black wheels PU 10' D250x40 x 10' D250x40

10. Transport

The rollator is a foldable product, convenient for transport. It is recommended to transport it in the folded position, as described above. During transport, ensure the product is in a stable position and cannot move around.

11. Service

- Always check that all parts are in place and securely fixed.
- Make sure all screws are always tightened.
- Immediately replace any worn or malfunctioning parts.

12. Reuse

The rollator is suitable for reuse. Before being passed on to another user, it must undergo a technical inspection at an authorised service centre.

Inspection points:

- Visual inspection of the frame
- Secure fastening of nuts and bolts
- Condition of plastic parts
- Condition of wheels and ball bearings
- Effectiveness and condition of the brakes

If any part is damaged or broken, it must be replaced with a new one. It is recommended to replace the wheels in every case of user change. All loose screws and nuts must be tightened. Before reuse, the product must be thoroughly cleaned and disinfected.

Warning! The rollator must not be reused if the frame is bent or broken. The user manual must always be provided with the product when changing users.

13. Disposal

At the end of its service life, the device should be handed over for recycling to a company specialising in disposal.

14. Warranty

The warranty does not cover improper use of the product. It excludes damage resulting from incorrect operation, accidents, and/or normal wear and tear. Any modifications to the product will void the warranty. Mechanical damage is not covered by this warranty.

For warranty repairs, please contact your local dealer.

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