PEDAL TRAINER
Art. 272002 & 272003

Gebrauchsanleitung
Brugervejledning
User Manual
Manuale
Bruksanvisning
Instrukcja Obsługi
Manuel de l'utilisateur

Pedal trainer **Art. 272002**

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Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

1. Introduction

Thank you for buying a high-quality Mobilex product. These instructions of use contain a description and important information about the use of the medical device. It is intended to help you to use the product correctly and safely. It is very important to read the instructions carefully before using the device. Pay particular attention to the safety guidelines and always follow them. Mobilex is constantly trying to improve and upgrade its products. Therefore we reserve the right to change specifications and technical data of the product without further notice.

Intended use of the device

The pedal trainer is a medical device which has been designed to be used for rehabilitation exercises that can be carried out by people with paralysis of the limbs. Thanks to its universal design, it enables both upper and lower limb exercises. Any other use is prohibited.

Indications

Due to the wide range of possible applications, indications cannot be directly linked to a specific clinical picture. Rehabilitation pedal trainers are indicated for use, as a rehabilitation exercise device, in the process of treating a number of diseases associated with severe limitations of movement of upper and / or lower extremities and a decrease in muscle strength. For the above reasons, it is necessary that the selection of the right rehabilitation device, suitable for a given clinical picture, should always be made by a specialist doctor or physiotherapist.

Contraindications

Contraindications cannot be related to one specific clinical picture and they depend on individual evaluation of potential user's health status. The choice of orthopedic aid appropriate for a given clinical picture must be always made with the help of a specialist physician or physiotherapist. In addition, observe the safety information in Chapter 2.

Areas of use

The product is suitable for indoor use and can be used by private individuals in domestic environments, in home care, in health and care facilities and in hospitals.

Medical incidents and unwanted events

In case of any unwanted events related to the device, please contact the dealer or manufacturer:

Mobilex A/S		
Grønlandsvej 5	Phone:	+45 87 93 22 20
DK-8660 Skanderborg	Telefax:	+45 87 93 17 77
Denmark	Email:	info@mobilex.dk

Safety requirements

- Make sure that this manual has been read by all persons using the device.
- The manufacturer assumes no responsibility for damage and/or injury resulting from failure to follow the operating instructions.
- Only use the product in a flawless condition.
- If defects or errors are detected you must immediately stop using the device and contact the dealer.
- Follow all manufacturer's recommendations and warnings.
- Be aware of the information on the product label.
- Only use the product for the described purpose.
- Don't stand on the pedal trainer
- Avoid making constructive changes to the device, unless you have the manufacturer's written acceptance for such modifications.
- The device must not be loaded with more weight than indicated in Technical data.
- The product should be completely assembled and unfolded before use. Please contact us if there are any problems in this regard.
- Folding and adjusting mechanisms can become a cause of danger if body parts get into the gap between the components and could become trapped when the gap narrows.
- Always check that the product is well unfolded before each use.
- Clean the product regularly.
- Do not use the product on slippery or unstable surfaces.
- Use only while seated in a chair or on bed with the pedal trainer placed securely on the ground or, in the case of arm training, placed securely on a surface at a height suitable for the user.

Parts supplied

The device is supplied in a carton box. The contents of the package are:

- ✔ 1 Pedal trainer;
- ✔ 1 User manual;

Preliminary check

Check the contents of the box with the list above. If at this stage you discover that a part is missing or damaged, contact your supplier immediately.

Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

3. Description

Product description

The Pedal trainer is made of 2-coloured steel frame. The pedals have safety straps and the pedal resistance can be adjusted. Due to its small size it can be set up anywhere, and thanks to its easy folding function it is easy to use. Appropriately shaped anti-slip caps on the base feet prevent the pedal trainer from sliding during exercise.

The Pedal trainer consists of the following elements:

1. Two pedals
2. Two foldable legs
3. Anti-slip caps
4. Knob for adjusting resistance
5. LCD display (only **Art. 272003**)
6. Steel frame
7. Locking button for legs in unfolded position

Assembly description

This device doesn't require an assembly instruction.

Function description

Using the pedal trainer

Spread the legs of the pedal trainer outward so that, thanks to the locking buttons, they snap into the unfolded position. The leg is unfolded when you hear a distinct click and it cannot be folded back. To fold the legs, first press and hold the locking buttons. The pedal trainer can be placed on the floor, bed, table, etc., depending on your needs. Using the knob, you can smoothly adjust the amount of resistance you need to overcome during exercise. The size of the safety straps on the pedals can be increased or decreased depending on the size of your feet/hands thanks to the adjustment fasteners.

WARNING! The pedal trainer may get hot during prolonged use!

Use of the optional LCD display (Art. 272003)

The LCD display shows the training parameters.
Button: **MODE/RESET** – a short press allows you to switch the display between available functions. Holding the button for 4 seconds resets the measured parameters.

- **TIME**: determines the time the pedal trainer has been operating.
- **CNT**: determines the number of rotations performed since the display was reset to zero.
- **CAL**: determines the estimated number of calories burned during exercise.
- **RPM**: determines the exercise speed (in rotations per minute).
- **SCAN**: starts the cycle of automatic presentation of all measured values in sequence.

Battery replacement*:* if the display stops working, replace the battery. To do this, remove the display from the pedal trainer structure and remove the cover from the bottom of the display. Replace the battery, close the lid and attach the display to the pedal trainer in the same place where it was before.


Accessories

There are no accessories provided for this device.

4. Technical data

Pedal trainer	Art. 272002	Art. 272003
Dimensions (unfolded)	44 x 39,5 x 23 cm	44 x 39,5 x 23 cm
Dimensions (folded)	44 x 27 x 14 cm	44 x 27 x 14 cm
Weight	2,16 kg	2,45 kg
Max. user weight	100 kg	100 kg

Product label

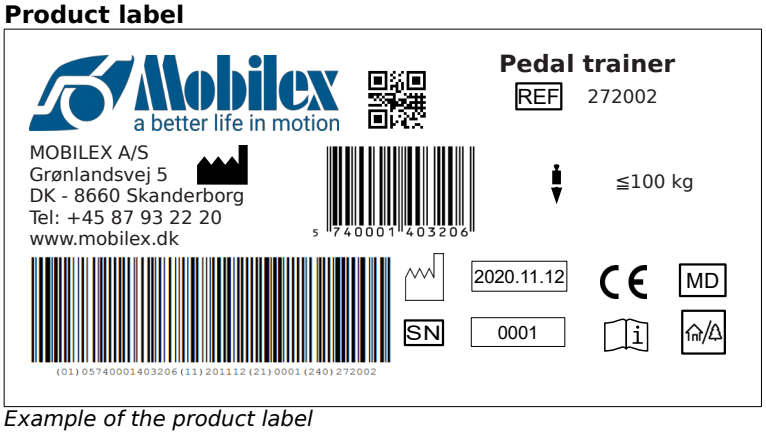Example of the product label

	QR-code link to homepage		Manufacturer
CE	CE marked product	REF	Item Number
	Observe the user manual		Max load
MD	Medical Device		Manufacturing date
	To be used indoor and outdoor	SN	Serial no

Meaning of the icons

Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

5. Care and maintenance

Regular dirt can be removed with standard cleaning products and a sponge or soft rag. Check the specific product information and use only commercial cleaning products which are suitable for cleaning. For disinfection, use only disinfectants certified for use with medical devices available locally. Do not use solvents or abrasives.

Despite the solid construction and the use of resistant materials the product is subject to wear. It is therefore recommended to check the product at regular intervals by a professional service. We recommend the following:

- ✔ Cleaning and disinfection of the product
- ✔ Checking the frame for bends, damage, wear or corrosion
- ✔ Checking wear parts and replacing them if necessary
- ✔ General re-tightening of the components

Recommended service life: Normally 10 years, depending on use.

Storage: We recommend storage in a dry place in temperatures above 0 °C.

Transportation

It is recommended to fold the device for transportation. The device should be transported in a stabilized and secured way with no risk of sudden or dangerous movements during transport.

Reuse

The product is suitable for reuse. Before a new user takes over the product, the device must be checked on a technical level and disinfected according to point 5.

6. Most common problems and solutions

Problems	*Possible causes*	*Solutions*
Pedal trainer isn't rotating	Resistance is set too high	Lower the resistance using the adjustment knob
Pedal trainer is wobbly	Pedal trainer isn't standing evenly on 4 legs	Put the pedal trainer evenly on a flat surface
LCD screen doesn't work	Battery is dead	Replace the battery

7. Disposal of the product

The product can not be disposed of with household waste and must be brought to the local recycling center.

8. Guarantee terms

Mobilex A/S offers 2 years of warranty for damages caused by errors in production or material errors. The warranty is not valid by not-intended repair or usage of the product. Parts subject to normal wear are not covered by the warranty, unless wear is caused by a manufacturing fault. If product changes are made without our written acceptance the CE-marking and product warranty will not be valid. In the event of product damages covered by this warranty, please notify your dealer or Mobilex A/S directly. The warranty does not cover transport costs and does not include compensation for personal injuries or non-use of the product during repair. The warranty does not cover damages due to the user manual not being followed.

9. Declaration of conformity

The device has passed the conformity assessment in accordance with EU Regulation: MDR 2017/745 regarding medical devices and is CE marked. A declaration of conformity has been issued for the device and is available at mobilex.dk/medias/2021-07/ce_pedaltraener__272002_3__8254.pdf or via the QR code on the right.

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