Bedienungsanleitung
Instruction manual
Manual de instrucciones
Mode d´emploi

S-VR

DE

EN

ES

FR

Enjoy mobility.

DE

Inhalt

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EN

Table of content

1. Preface 40
1.1. Signs and symbols 40
2. Important safety instructions 41
2.1. General safety instructions 41
2.2. Risk of tipping over 43
2.3. Safety in vehicles 43
2.4. Use on Public Roads 43
2.5. Brakes 43
2.6. Purpose 44
2.7. Indications 44
2.8. Contraindications 45
2.9. Declaration of conformity 45
2.10. Liability 45
2.11. Service life 45
3. Description of Product and Con- 46
3.1. Inspection on Delivery 46
3.2. Contents 46
4. Type plate and serial number 47
5. Summary 48
6. Assembly/Adjustments 49
6.1. Folding and unfolding 49
6.2. Legrests 50
6.2.1. Fitting/removing legrests 50
6.2.2. Height adjustment of the calf 51
6.2.3. Angle adjustment of the leg- 51
6.3. Adjusting the legrest 51
6.4. Sideguards 52
6.4.1. Removing the sideguards 52
6.4.2. Fitting the sideguards 52
6.5. Backrest 52
6.5.1. Angle adjustment of the backrest 53
6.6. Backrest extension/ 53
6.7. Wheels and tyres 53
6.8. Brakes 54
6.8.1. Operating the parking brake 54
6.8.2. Adjusting the parking brake 55
6.8.3. Drum brake (optional) 56
6.8.4. Anti-tip wheels and tipping aids 56
7. Driving 57
7.1. Getting in and out of the chair 57
7.2. Propulsion 58
7.3. Going up or down stairs or a 59
7.4. Driving the wheelchair 60
8. Accessories 60
8.1. Lap belt (optional) 60
8.2. Therapy table (optional) 61
8.3. Infusion/oxygen bottle 61
9. Technical Data 62
10. Transport 64
10.1. Transporting the wheelchair 64
10.2. Transportation in 64
11. Care instructions for the user 65
11.1. Cleaning and disinfection 65
11.2. Check Before Use 65
11.2.1. Checking the brakes 65
11.2.2. Checking the tyres 65
11.2.3. Oiling the moving 66
12. Repairs 66
13. Maintenance instructions 66
14. Use by others 66
15. Storage/Shipping 67
16. Disposal 67
17. Warranty 68

1. Preface

Dear user,
You have decided on a high-quality Bischoff & Bischoff wheelchair, and we would like to thank you for your trust.
This Instruction manual contains all the information necessary to operate the S-VR wheelchair.
Your wheelchair is designed for both indoor and outdoor use.
Please read and follow the Instruction manual carefully before using your new wheelchair for the first time. The Instruction manual is an integral and necessary component of the wheelchair. Keep the Instruction manual where you can refer to it and if you transfer the wheelchair, give the manual to the new owner.
For visually impaired users, this document is available as a PDF file at our website www.bischoff-bischoff.com.

Note!
In spite of careful research and editing, incorrect information may have crept into our Instruction manuals. We reserve the right to make technical modifications. Illustrations may include accessories. The current version of each Instruction manual is ready for download on our website:
www.bischoff-bischoff.com

1.1. Signs and symbols

Important! This is especially useful information on a given subject.

Warning! These are especially safety-relevant Instructions. Follow the directions in the Instruction manual!

2. Important safety instructions

2.1. General safety instructions

- To avoid falls and dangerous situations, first practice driving your new wheelchair on level, unobstructed ground. An attendant is recommended.
- Before operating the wheelchair, check that all attached parts are correctly fastened.
- Make sure that the screw axles are located properly in the axle holders in order to prevent a risk of injury.
- Only lift the wheelchair using firmly installed parts (see Chapter 7.3).
- The risk of tipping over to any side increases with changes in balance through body movements or loading the wheelchair.
- The wheelchair may not be used for transport of more than one person or any additional load. Please respect the purpose of the product (see Chapter 2.6)!
- When carrying out setting or adjustment work on the wheelchair there is a risk that body parts could get pinched!
- Keep items of clothing and body parts away from gaps between moving parts!
- When going up or down obstacles (steps, etc.) call in an attendant.
- Only go up or down steps or large obstacles with the help of at least two attendants (see Chapter 7.3). Use equipment such as ramps, lifts or stair-climbing devices, if these are available. When going up or down obstacles, the anti-tip wheels can be used as tipping aids.
- When driving on downhill/uphill slopes, please be aware that even at a few degrees, there is a risk of tipping over. We recommend the use of anti-tip wheels (see Chapter 6.8.4).
- Avoid collisions of any type! Avoid hitting obstacles (steps, kerbs) at full speed or jumping down steps.
- Only drive on a level, firm surface.
- When transferring to or from the wheelchair, apply both parking brakes (see Chapter 6.8.1).
- Replace the tyres as soon as the tread is worn down.

- When on roads, observe the road traffic regulations (StVZO).
- Each time, before setting off, check the reflectors of your wheelchair to make sure they are not damaged and are visible for other road users.
- When using the handrims as operational brakes (braking from speed or on long downhill slopes) your fingers can quickly heat up. We recommend using leather gloves when driving outdoors, as these will increase the grip and protect your fingers at the same time.
- Only use the parking brakes for parking, they must not be used to brake when you are driving.
- The seating materials are flame resistant. However, please keep cigarettes and other possible ignition sources away from your product.


Protect the wheelchair from sunlight to avoid burns caused by heated components

- Observe the storage and operating conditions (see Chapter 15).
- When transferring between the wheelchair and the bed (see Chapter 7)
 - apply both parking brakes
 - remove the relevant sideguard (see Chapter 6.4)
 - Flip the footplates up (when transferring to or from the chair) or remove the legrests (see Chapter 6.2)
- Reporting incidents. The operator or user must immediately report
 1. any malfunction
 2. any change to the features or performance as well as
 3. Immediately report any improper labelling or instructions for use of a medical device that have or could have led to the death or serious deterioration of the health of a patient, an employee or a third party to your competent regional authority for medical devices. This institute will immediately forward the notification to the operator's competent authority and also inform the manufacturer and the manufacturer's competent authority.

EN

2.2. Risk of tipping over

The anti-tip wheels which are fitted on both sides (see Chapter 6.8.4) will largely prevent the wheelchair from tipping over backwards. However please note that driving with anti-tip wheels is only possible in a limited way on ground with different levels, lift platforms and ramps. Make sure there is sufficient gap at the top, bottom and sides.

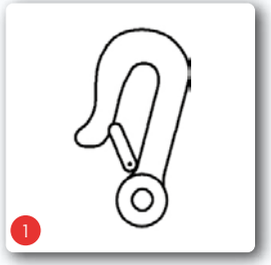2.3. Safety in vehicles

To see whether your wheelchair is permitted for use as a seat in a vehicle for transporting individuals with impaired mobility or not, you can refer to the sticker next to the type plate on the wheelchair.

- Whether your wheelchair is permitted for use as a seat in a vehicle for transporting individuals with impaired mobility or not (Fig. 1)
- Whether your wheelchair is not permitted for use as a seat in a vehicle for transporting individuals with impaired mobility or not (Fig. 2)

The S-VR wheelchair is not suitable for use as a car seat,
so it may not be used as such, under any circumstances. It does not meet the requirements stipulated by ISO 7176-19. Failure to observe this instruction can lead to extremely serious or fatal injuries in the event of an accident.

2.4. Use on Public Roads

The wheelchair is designed for indoor and outdoor use. Please keep in mind that you are using public roads and therefore have to observe all road traffic regulations (StVZO). Do not endanger other road users through inconsiderate driving, especially on walkways.

2.5. Brakes

Operating brake:
When in motion, brake the wheelchair using the handrims. Please note that the hands may become hot when doing this.

(see Chapter 6.8)
Parking brake/parking:
As soon as the wheelchair has stopped, push both brake levers forward. (see Chapter 6.8.1). The wheelchair is securely braked when it can no longer be pushed and both brake levers are pushed forward to the stop position.
Pull the lever towards you to release the brake (see Chapter 6.8.1).
Drum brake:
As an option, the wheelchair can be fitted with a drum brake which can only be operated by an attendant using the brake levers on the push handles

2.6. Purpose

The wheelchair is designed to increase the mobility only for adults who have difficulty walking with the indications specified (see Chapter 2.7). The maximum user weight of 120 kg for the S-VR must not be exceeded.
Please note:
We can only provide a warranty if the product is used under the conditions specified and for the intended purposes.

2.7. Indications

Provision of a wheelchair is intended for adults who are unable to walk or who have severely impaired mobility, e.g. through:

- paralysis
- loss of limb
- defect/deformity of limb
- joint contractions (not on both arms)
- joint damage (not in both arms)
- other illnesses.
 Provision of reinforced wheelchairs with backrest adjustment up to 30 degrees is advised, if the body weight of the wheelchair user exceeds 90 kg, normal wheelchairs no longer offer the required safety and because of the disability, an almost horizontal position must be adopted or is necessary from time to time.

2.8. Contraindications

The provision of a wheelchair is unsuitable for persons with:

- no capability to sit
- joint damage/contractions on both arms
- perceptual disorders
- severe equilibrium disorders
- loss of limbs on both arms
- impaired or insufficient vision
- strongly impaired cognitive abilities

2.9. Declaration of conformity

We Bischoff & Bischoff GmbH declare under our sole responsibility, that the S-VR wheelchair meets all the requirements of MDR 2017/745 that are applicable.

2.10. Liability

We can only provide a warranty if

- the product is used under the conditions specified and for the intended purpose.
- Modifications, enhancements, repairs and maintenance work is carried out only by those individuals who have been authorised by us to do so, and
- the wheelchair is used in accordance with the instructions for use.

2.11. Service life

The expected product service life is up to 5 years. Prerequisites for this are proper use as well as compliance with the safety, care and maintenance information in accordance with this Instruction manual.

3. Description of Product and Contents

3.1. Inspection on Delivery

All Bischoff & Bischoff GmbH products undergo a proper final inspection in our company and are labelled with the CE mark (see Chapter 4).
The wheelchair will be delivered in a special cardboard box. Keep the box after unpacking, if possible. It is useful for later storage of the product or for sending it back.
Delivery and first instructions are taken care of by your qualified medical supply dealer.
If the scooter is shipped by rail or road, the goods should be checked immediately for transport damages in the presence of the bearer.
Check the contents for completeness (see Chapter 3.2) and make sure nothing is damaged. Should you find any irregularities or damage, please contact our customer service.

3.2. Contents

After receipt of the goods, please check the contents for completeness right away. The contents consist of:


- Outer packaging.
- wheelchair (preassembled),
- legrests
- this Instruction manual
- any accessories

EN

4. Type plate and serial number

The type plate (Fig. 3) and serial number are located on the cross brace.

- A Manufacturer logo
- B Manufacturer information
- C MD Medical device
- D CE mark
- E Observe safety in vehicles
- F Warning! - read the User manual
- G Important! - read the User manual
- H Model designation
- I Model number
- J Maximum user weight
- K Date of manufacture
- L UDI number
- M Serial number

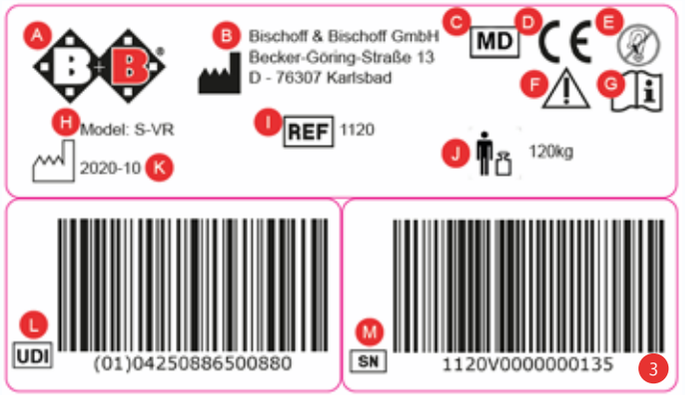Depending on the model number, the information may vary. Refer to your product's type plate for the information which relates to your product.

5. Summary

(see Fig. 4)

- A Push handle
- B Back cushion
- C Front wheel
- D Sideguard
- E Seat pad
- F Drive wheel
- G Quick-release axle
- H Parking brake
- I Handrims
- J Legrests
- K Backrest extension (S-VR)
- L Backrest adjustment (infinite)
- M Armrest pad
- N Calf pad
- O Footplate (flip-up)

4

EN

6. Assembly/Adjustments


Your new Bischoff & Bischoff wheelchair will be delivered in a Bischoff & Bischoff original cardboard box, assembled and folded.
In order to avoid transport damage, detachable attachments are delivered separately with the product.
Please see Chapter 5 for information on which accessories are available for your wheelchair.

6

6.1. Folding and unfolding

Unfolding


- Position the folding wheelchair next to you.
- Tip the wheelchair towards you until the weight is off one of the drive wheels.
- Push the seat tubes apart and push them down with the flat of your hand, until they lock into place in the bracket on the frame (Fig. 5).
- The seat pad must be completely unfolded.
- Apply the left and right parking brakes on the drive wheels.
- Fit the legrests and flip them inwards (see Chapter 6.2). Flip the footrests down (Fig. 5)
- Swivel the cross tube upwards and fix it on both sides using the star screws (Fig. 6 A/B).
- Fit the backrest and fix it on both sides using the star screws (Fig. 6 C)
- Now you can sit down

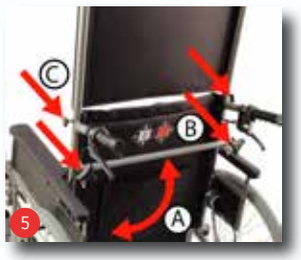5

Folding

- First remove the backrest extension by undoing the star screws on both sides and pulling the backrest extension upwards (Fig. 7 A).
- Undo the star screws of the cross tube on both sides and swivel it downwards (Fig. 7 B/C).
- Flip the footplates and calf pads upwards (Fig. 8).
- Remove the seat cushion, take hold of the front and rear of the seat pad and pull it up,

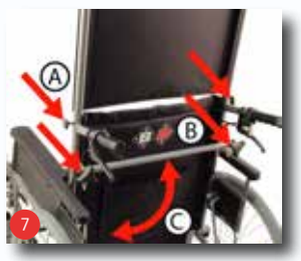7

8

When folding and unfolding, watch your fingers as there is a risk of pinching

6.2. Legrests

6.2.1. Fitting/removing legrests

You can flip the footplates up separately to the side and adjust to the lower leg length of the user (see Chapter 6.3). You can remove the legrests, swing them away and refit them without needing to use tools.

Attach legrests

- The legrests are inserted into the receivers on the frame (Fig. 9).
- Then swivelled inwards with a gentle push. Done.
- To disassemble the legrest, pull the lever (Fig. 10) on the inner side, swivel the legrest back out and pull it out of the receiver.

Legrest with footplate

Both legrests are delivered with footplates fitted.
When seated, the footplates on both legrests can be flipped up (Fig. 8).
Before setting off, fit the legrests as described above.

Never stand on the footplate and when swinging it away, watch your fingers as there is a risk of pinching

EN

6.2.2. Height adjustment of the calf pad

- Flip the calf pad forwards (Fig. 11).
- Move it to the height required. 3 heights are possible
- The calf pad is then swivelled back.

6.2.3. Angle adjustment of the leg-rest

The legrests can be adjusted separately up to 90°. Push the locking lever down (Fig. 12).

- swivel the legrest by simply pulling it upwards to the required position and then let go of the lever
- the legrest will lock into place automatically
- To move the legrest into a lower position or the initial position, pull the locking lever and push the legrest back to the required position.

To change seat, the legrest must be swung away or removed. Never sit on the legrest. Risk of tipping over!

6.3. Adjusting the legrest length

The length of the legrest can be infinitely adjusted using the telescopic tube.

- Undo the hexagon head screw (SW13) by turning it to the left (Fig. 13).
- Adjust the legrest as required
- Re-tighten the screw by turning it to the right.

To change seat, the legrest must be swung away or removed. Never sit on the legrest. Risk of tipping over!

6.4. Sideguards

The sideguards on the S-VR are detachable.

6.4.1. Removing the sideguards

- Rotate the star screw on the handlebar stem of the wheelchair to the left to undo it (Fig. 14).
- Now pull the sideguard upwards to remove it (Fig. 14).

6.4.2. Fitting the sideguards

- Insert the sideguard into the frame receiver at the front and back (Fig. 15).
- Re-tighten the star screw by turning it to the right.

When fitting the sideguards, make sure that clothing or other items cannot be trapped. Before using for the first time, check that the sideguards are firmly in place!

Never use the sideguards to carry the wheelchair!

6.5. Backrest

The angle of the backrest can be infinitely adjusted from -4° forwards to 30° or 82° backwards. The conversion from 30° to 82° or vice versa is not carried out at the factory, as the HMV no. is only valid for backrest angles up to 30°.

EN

6.5.1. Angle adjustment of the backrest

- Apply the parking brakes
- Press both Bowden cable levers of the compressed spring system and keep them pressed (Fig. 16)
- Adjust the backrest angle as required
- Release the Bowden cable lever (Fig. 16).

6.6. Backrest extension/ Headrests

As standard, the S-VR wheelchair is fitted with a large backrest extension which can be used as a headrest. This is inserted from above into the guide provided and secured with the star screw (Fig. 17).

6.7. Wheels and tyres

Drive wheels

PU tyres are fitted on your S-VR as standard. When worn or the tread is worn down or defective, please ask your specialist dealer to service them.

Swivel wheels

PU tyres are fitted on your S-VR as standard. When worn or the tread is worn down or defective, please ask your specialist dealer to service them.

If the drive wheel or front wheels are worn or defective, ask your dealer to service them!

6.8. Brakes

As standard, the wheelchair is fitted with parking brakes to be operated by the wheelchair user.

The parking brakes must not be used as operational brakes, as there is a risk of death due to the risk of falling!

Please note that there is a risk of burns if your hands heat up when braking using the handrims!

6.8.1. Operating the parking brake

The parking brake must only be used for parking. To slow down when driving you must only use the drive wheel handrims.


- When braking using the handrims, hold the highest point. First let the handrims slide through your hands.
- Slowly increase the braking force by exerting gripping force on the handrims, until you come to a standstill.
- As soon as the wheelchair has stopped, push both brake levers forward (Fig. 18). The wheelchair is securely braked when it can no longer be pushed and both brake levers are pushed forward to the stop position.
- Pull the brake lever backwards again to release the brakes.

EN

6.8.2. Adjusting the parking brake

Adjusting the brakes must only be carried out by trained personnel. Please bear in mind that if the brakes are set incorrectly this can have life-threatening consequences.

- Parking brakes only work on one drive wheel at a time! Adjustments must therefore be made on both wheels separately!
- If you should notice that the brakes are not working properly or evenly, proceed as follows:
- Check the tyre tread depth. Also make sure that the wear is even. If the tread depth in the centre is worn down, the tyres must be replaced by a specialist dealer.
- Check the brake setting. This must be adjusted so that the gap between the tread of the wheel and the brake pins is 16 mm when the brake is open (Fig. 19 / A).
- Once you have finished the adjustment, tighten the screws you have loosened and carry out a brake test!

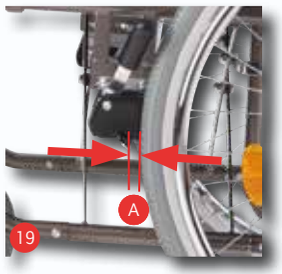Adjustments to the brake system should only be carried out by an authorised dealer.

Never use the brake lever as a support - e.g. when transferring or standing up -, otherwise there is a risk of tipping over and falling! The lever could break off!

The parking brake - with or without brake lever extension - only works on one drive wheel at a time and must therefore always be used on both sides!


6.8.3. Drum brake (optional)

As an option, the wheelchair can be fitted with a drum brake which can only be operated by an attendant using the brake lever on the push handles (Fig. 20).

Adjusting the drum brake

Basic setup for the drum brake is carried out on the cable clamping screw (A). Fine adjustment is carried out using the adjusting nut (B) on the brake cable. Unscrewing the adjustment nut will increase the braking effect (Fig. 21).

- Drum brakes only work on one drive wheel at a time!
- Drum brakes are also subject to natural wear. If the braking effect decreases, the wheelchair must no longer be used. In order to avoid this happening, for your own safety you should check that the brakes are working properly each time before setting off!
- Once you have finished the adjustment, re-tighten all screws and carry out a brake test!

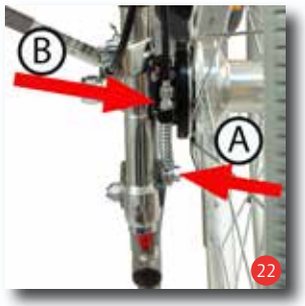Adjustments to the brake system should only be carried out by an authorised dealer.

6.8.4. Anti-tip wheels and tipping aids

To increase tipping stability on each side there is a tube with a small wheel. They prevent the wheelchair from tipping over backwards. The anti-tip wheels can also be used as a tipping aid by the attendant.

Activate anti-tip wheels

- Press the locking button in the tube (Fig. 22)
- Rotate the tube such that the wheel is pointing downwards in the direction of the road. The tube can be locked in two extension lengths.
- If you do not need anti-tip wheels and would like to use these as tipping aids, proceed as described above, just then rotate the wheels so that they are pointing upwards (Fig. 23).

"When" going up or down obstacles (e.g. kerbs), the anti-tip wheels must be rotated by 180°
.

With the wheels pointing upwards, the tipping protection is disabled; an attendant is recommended (Fig. 23).

For safety reasons, always use the anti-tip wheels as a pair!

7. Driving

7.1. Getting in and out of the chair

Side transfer:

Never stand on the footplate and when swinging it away, there is a Risk of tipping over!

- Move the wheelchair as close as possible to the seat surface you are sitting on.
- Secure the wheelchair by applying the parking brake.
- If you are transferring from another wheelchair or indoor frame, secure these by applying the parking brakes.
- Now remove the sideguard and legrest and slide sideways onto the seat surface of the wheelchair. Please make sure that you sit as far back in the seat as possible.
- Then refit the armrests and legrests.

Frontal transfer:

Never stand on the footplate and when swinging it away, there is a risk of tipping over!

- Move the wheelchair as close as possible to the front of the seat.
- Secure the wheelchair by applying the parking brake.
- If you are transferring from another wheelchair or indoor frame, secure these by applying the parking brakes.
- Flip the footrests up.
- Now slide onto the seat surface by rotating your body. Please make sure that you sit as far back in the seat as possible.
- Then flip the footrests back down and put your feet on them.

Before transferring on or off the wheelchair apply the parking brake. You must only drive outside your home if you have sufficient driving experience. Where possible, you should have an attendant.

7.2. Propulsion

Propulsion is achieved using the handrims on the drive wheels, as long as it is not being pushed by an attendant. Carry out the first test drives carefully, until you have got used to the wheelchair and its driving characteristics. Wheelchairs only have limited security against tipping over and slipping.

Take extra care when braking, starting off or turning on slopes. Extended leaning out of the wheelchair also increases the risk of tipping over. For greater safety, all wheelchairs are fitted with a lap belt (not standard).

In order to be able to go up or down obstacles more easily, the attendant has the option of lining the wheelchair up at right angles to the obstacle by pressing their foot on the tipping aid.

In accordance with StVZO 66/1 of 01.01.1981, manual wheelchairs which are used on public roads (on pavements too) require lights which

are permanently fitted.

This consists of (passive):

- two red reflectors at the back.
- Reflectors in the spokes of the drive wheels.

If the wheelchair is used after dark, then in addition, lights are required (active):

The functionality of lights must be guaranteed by a circuit which is independent of the drive system (batteries).

7.3. Going up or down stairs or a high step

Only the frame and the push handles must be used as grip points.

The helpers only grip on the firmly fixed frame parts (Fig. 24). In order to go up stairs, the rear helper pulls the wheelchair backwards up the steps by holding onto firmly-mounted handles. The lower helper grips on the front frame tubes and stabilises the wheelchair. When doing this he pushes the drive wheels of the wheelchair onto the steps. When descending, the lower helper brakes the chair by pressing the wheelchair straight into the steps. The upper helper holds the wheelchair by the push handles, keeps it secure and holds it in the right position. The wheelchair should roll down step by step.

Two attendants must help with obstacles and stairs with more than three steps.

7.4. Driving the wheelchair

- Sit in your Bischoff & Bischoff wheelchair. Make sure that you sit as far back as possible.
- Propulsion is achieved using the handrims on the drive wheels. Place your hands on the handrims using your thumbs and bent index fingers.
- From this position, push the wheelchair with both hands.
- Turning right: hold the right handrim tight and only push with the left one
- Turning left: push in the opposite way
- Turning on the spot: push the handrims in opposite directions at the same time
- Braking: reduce the speed of the handrims with your hands.

Carry out the first test drives carefully, until you have got used to the wheelchair and its driving characteristics. Wheelchairs only have limited security against tipping over and slipping. Take extra care when braking, starting off or turning on slopes. Extended leaning out of the wheelchair also increases the risk of tipping over. Tipping protection reduces this risk. In order to go up or down obstacles more easily, the attendant can tip the wheelchair (see Chapter 6.8.4).

Please note that there is a risk of burns if your hands heat up when braking using the handrims!

8. Accessories

8.1. Lap belt (optional)

You can order a lap belt as an option and this will secure users who do not have the required stability in the wheelchair.
This is fitted to the frame of the wheelchair backrest with the sideguard temporarily removed.

- Wrap the open end of the belt once round the frame tube of the back section and then fix it with the buckle
- Guide this through the backrest on both sides. You can also adjust the length of the belt using the buckle. Adjust the belt length so that the belt is taut but does not pinch
- Use the fastening clip to open and close the belt.

8.2. Therapy table (optional)

The table is simply pushed into the armrest from the front. To remove the therapy table, simply pull it back out.

8.3. Infusion/oxygen bottle holder (optional)

The infusion/oxygen bottle holder is fixed on the backrest tube of your wheelchair.
Warning - with increased load from the oxygen bottles, the risk of your wheelchair tipping over backwards will also increase. We strongly recommend the use of anti-tip wheels.

Please note that reflectors must not be covered.

Warning - with increased load in the bag, the risk of your wheelchair tipping over backwards will also increase. We strongly recommend the use of anti-tip wheels.

EN

9. Technical Data

	S-VR
Dimensions and weights	
Total length with legrests (mm):	1265
Total length without legrests (mm):	940
Overall width (mm):	700
Total height (mm):	880
Total height with backrest extension	1330
Kerb weight (kg):	25.4
User weight (kg):	120
Default seat angle (°):	4°
Backrest angle	- 4° to + 30° (without HMV approval up to + 82°)
Seat width (mm):	400, 440, 480
Seat depth (mm):	420
Front seat height without seat cushion (mm):	500
Backrest height without seat cushion (mm):	430
Backrest height with backrest extension	840
Armrest height without seat cushion (universal sideguard) (mm):	220
Armrest length (mm):	545
Armrest length to backrest tube	410
Tyres	
Drive wheels (inch):	24" x 1 3/8" PU
Front wheels (mm):	200 x 50
Handrim diameter (mm)	530
Handrim tube diameter (mm)	20
Wheelbase (mm)	515

	S-VR
Camber (°)	0°
Material	
Frame:	Steel
Cross brace, sideguard tube, legrest tube	Steel
Seat and backrest:	Dartex
Storage conditions	
Temperature:	-40° - +65° C
Air humidity	45 - 60 %
Operating conditions	
Temperature:	-25° - +50° C
Air humidity	20 - 80%
Accessories and options for all models	
Lap belt, therapy table, infusion/oxygen bottle holder.	
Accessories and options for individual models	

The dimensions stated are below the manufacturing tolerances of +/- 10 mm.

10. Transport

10.1. Transporting the wheelchair

In order to make the wheelchair as manageable as possible to transport, remove all detachable parts (legrests, drive wheels, armrests, loose accessories) and fold the wheelchair up
(see Chapter 6.1).

10.2. Transportation in vehicles

The S-VR wheelchair is not suitable for use as a car seat. Under no circumstances may it be used as such. It does not meet the requirements stipulated by ISO 7176-19. Failure to observe this instruction can lead to extremely serious or fatal injuries in the event of an accident.
The symbol (Fig. 25) on the type plate on the wheelchair indicates that it must not be used as a vehicle seat.

If your wheelchair is not crash-tested, under no circumstances must it be used as a vehicle seat in a vehicle for transporting individuals with impaired mobility. There is danger of fatal injury!

EN

11. Care instructions for the user

11.1. Cleaning and disinfection

- Wipe the frame parts of the wheelchair with a damp cloth. If there is heavy contamination, also use a mild cleaner.
- Clean the wheels with a damp brush with plastic bristles (do not use a wire brush!).
- Wash the push handles, armpad, seat and back cushion with a mild cleaner.
- For disinfection use a disinfectant listed by the Association for Applied Hygiene.
- Upholstery parts may be machine washed at 40°C. Do not machine dry, dry at room temperature only.

11.2. Check Before Use

11.2.1. Checking the brakes

Make sure the brake system is functioning properly before each use.

- With the parking brakes applied, the wheelchair cannot be pushed on dry surfaces. When the brakes are released, the wheelchair runs without scraping noises and good directional stability. If this is not the case, notify your dealer immediately and stop using the wheelchair until it has been serviced!
- Brake lever is firmly in place and cannot be twisted.
- All screws are tight.

The wheelchair must be checked in accordance with the maintenance plan
and be in proper technical condition.

11.2.2. Checking the tyres

Each time before setting off, check the tyres for damage and if you have pneumatic tyres check that the air pressure is correct.

11.2.3. Oiling the moving components

- All pivot points of the cross brace
- Moving parts of the sideguards
- Brake lever bearings
- Legrest bearings

Carried out by the user or an assistant. Before oiling the components remove any residual oil. Remove excess oil.

- Check that the screws are tight.

12. Repairs

Repair and servicing work must only be carried out by authorised partners (medical suppliers or dealers).

13. Maintenance instructions

Typical wear parts are brakes and tyres. Contact your dealer for spare parts.

We recommend that you have your wheelchair inspected by your dealer at least once a year.

If there are faults or defects, take the wheelchair to your dealer immediately for repair.

The information and documentation required for repair and maintenance shall be supplied to the mobility dealer on request. If you have questions, you can contact the Bischoff & Bischoff customer service.

14. Use by others

The wheelchair is designed for reusage by other persons. Before it is passed onto another user, it must be serviced by a dealer and hygienically cleaned. When passing on the wheelchair, don't forget to hand over this Instruction manual and the address of your dealer to the new user.

15. Storage/Shipping

If the wheelchair is to be stored or shipped, all loose parts which are not fixed, except for the drive wheels, must be removed and packed individually in an appropriate cardboard box. The individually packed parts can then be packaged together in a larger cardboard box. We recommend keeping the original packaging and storing it with these purposes in mind so that it is available if required. Your wheelchair then has optimum protection against environmental influences during storage or transport.

16. Disposal

If your wheelchair is no longer being used and is to be disposed of, contact your mobility dealer.

If you wish to dispose of the product yourself, contact local recycling companies to enquire about the regulations that apply to your place of residence.

17. Warranty

1. For products supplied by Bischoff & Bischoff, the period of warranty is 24 months from purchase date. The products are delivered free from production and material defects. In case of evidence of material or manufacturing defect, defective parts will be replaced free of charge.
2. This warranty excludes claims for the following reasons:
a) normal wear and tear, which include but is not limited to the following parts where fitted; Batteries, motor brushes, handles, armrests, upholstery, tyres, brakes, caps, etc.
b) overloading the product such as, for example, by exceeding the maximum user weight or load.
c) the product or part has not been maintained or serviced in accordance with the Instruction manual, the care and hygiene instructions or the manufacturer's recommendations as shown in the maintenance instructions.
d) Accessories have been used which are not original accessories.
e) The product or part has been damaged through negligence, accident or incorrect use.
f) changes/modifications have been made to the product or parts, which deviate from the manufacturer's specifications.
g) Repairs have been carried out, before our customer service has been informed of the circumstances.
3. To enforce the warranty, please contact your dealer with the exact details of the nature of the difficulty. The product must be repaired by a dealer authorised by Bischoff & Bischoff.
4. For parts, which have been repaired or exchanged within the scope of the warranty, the warranty will be extended to the remaining warranty period for the product in accordance with point 1).
5. For original spare parts which have been fitted at the customer's expense, these will have a 12 months guarantee, following the fitting.
6. If the repair fails after a reasonable time, the customer can, at his discretion, request a reduction, compensation or the termination of the purchase contract.
7. The guarantee is subject to the law of the country in which the product was purchased from Bischoff & Bischoff.

Please also observe the Bischoff & Bischoff care, warranty, hygiene and maintenance instructions. Your qualified mobility dealer will be happy to provide these to you.

EN

Warranty certificate

Product: S-VR / S-Care

Serial number:*

* (must be inscribed by the mobility dealer)

Dealer:

Date and seal

(Please note our general business terms and conditions)

*Signs are located on the cross brace

Note!

In spite of careful research and editing, incorrect information may have crept into our instruction manuals. We reserve the right to make technical modifications. Illustrations may include accessories. The current version of each Instruction manual is ready for download on our website:

www.bischoff-bischoff.com

Non-observance of the operating instructions as well as improper use, maintenance jobs and technical changes or adaptions without the permission of the Bischoff & Bischoff company, lead to guarantee loss and the loss of product liability in general.

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Bischoff & Bischoff France S.A.R.L.