SWIVEL SEAT CUSHION

Art. 278051, 278052, 278055, 278056

Gebrauchsanleitung
Brugervejledning
User Manual
Manuale
Bruksanvisning
Instrukcja Obsługi
Manuel de l'utilisateur

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Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

1. Introduction

Thank you for buying a high-quality Mobilex product. These instructions of use contain a description and important information about the use of the medical device. It is intended to help you to use the product correctly and safely. It is very important to read the instructions carefully before using the device. Pay particular attention to the safety guidelines and always follow them. Mobilex is constantly trying to improve and upgrade its products. Therefore we reserve the right to change specifications and technical data of the product without further notice.

Intended use of the device

The swivel seat cushion is a medical device designed to increase the mobility of persons who have difficulties moving and repositioning. The swivel seat cushion makes it easier for these people to get in or out of the car, bed or chair. The swivel seat cushion allows the user to sit on the swivel seat and easily turn using their own strength or be turned by an attendant. Any other use is prohibited.

Indications

Due to the wide range of possible applications, the indications cannot be tied to a specific clinical picture. The transfer products are indicated for multiple applications such as: positioning and/or transferring of patients with severely restricted mobility, positioning and/or transferring of patients who are highly sensitive to pain. Indications also depend on individual evaluation of potential user's health status by professional physician or physical therapist, but main target user groups are elderly people, persons with paresis, rheumatic diseases, and other people with limited mobility of the lumbar spine, knees or hips. Using the product for any purpose other than those described above is prohibited.

Contraindications

Contraindications cannot be related to one specific clinical picture, and they depend on individual evaluation of potential user's health status. The choice of orthopedic aid appropriate for a given clinical picture must be always made with the help of a specialist physician or physiotherapist. In addition, observe the safety information in Chapter **2**.

Areas of use

The product is suitable for indoor and outdoor use and can be used in home care, in health and care facilities and in hospitals.

Medical incidents and unwanted events

In case of any unwanted events related to the device, please contact the dealer or manufacturer:

Mobilex A/S		
Grønlandsvej 5	Phone:	+45 87 93 22 20
DK-8660 Skanderborg	Telefax:	+45 87 93 17 77
Denmark	Email:	info@mobilex.dk

2. Safety requirements

- Make sure that this manual has been read by all persons using the device.
- The manufacturer assumes no responsibility for damage and/or injury resulting from failure to follow the operating instructions.
- Only use the product in a flawless condition.
- If defects or errors are detected you must immediately stop using the device and contact the dealer.
- Follow all manufacturer's recommendations and warnings.
- Be aware of the information on the product label.
- Only use the product for the described purpose.
- Do not stand on the device.
- Avoid making constructive changes to the device, unless you have the manufacturer's written acceptance for such modifications.
- The device must not be loaded with more weight than indicated in Technical data.
- Always check that the product is in good condition without signs of tear and wear such as loose seams, frays and material inconsistency.
- Do not leave the product unattended on the floor, a person may accidentally step on it with the risk of injury.
- Always perform a risk assessment before starting the specific repositioning session in order to make sure the product model is correct and that the procedure is fully understood.

Parts supplied

The device is supplied in a carton box. The contents of the package are:

- ✔ 1 Swivel seat cushion;
- ✔ 1 User manual;

Preliminary check

Check the contents of the box with the list above. If at this stage you discover that a part is missing or damaged, contact your supplier immediately.

Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

3. Description

Product description
The swivel seat cushion adapts to most cars and surfaces and allows the user to sit on the swivel seat and easily turn using their own strength or be turned by an attendant.. The device is made of high quality polyester and PU memory foam, as well as a low friction material that provides good comfort and durability with silicone anti-slip technology.

Assembly description
This device doesn’t require an assembly instruction.

Function description


Using the swivel seat cushion
Place the cushion in the middle of the car seat (or other surface) with the dark, non-slip surface facing down. Stand with your back to the cushion. Bend your knees slightly by slowly moving backwards and sit on the cushion. Slowly turn your upper body towards the inside of the car and grab the steering wheel with your hands. Lift your legs together. Firmly hold the steering wheel and turn your legs and body into the car. Position yourself in the driver's seat and adjust the distance to the steering wheel. To get out of the car: hold the steering wheel, lift your legs and slowly turn your body out of the car by using the cushion. Put your feet on the ground and leave the car.

Accessories
There are no accessories provided for this device.

4. Technical data

Swivel seat cushion	Art. 278051	Art. 278052	Art. 278055	Art. 278056
Diameter ø	45 cm	50 cm	37 cm	44 cm
Max. user weight	100 kg	100 kg	130 kg	130 kg
Material	Polyester	Polyester	Polyester	Polyester


Product label

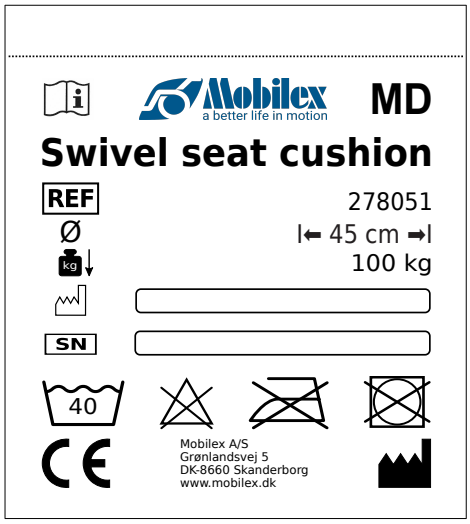Example of product label sewed on the swivel seat cushion

Meaning of the icons

	QR-code link to homepage	CE	CE marked product		Drip-dry
	Manufacturer		Manufacturing date		Medium ironing temperature
	Observe the user manual		Max load	30	Washing temperature
SN	Serial no.	MD	Medical Device		Do not tumble dry
REF	Item Number		Only to be used indoor		Handwash

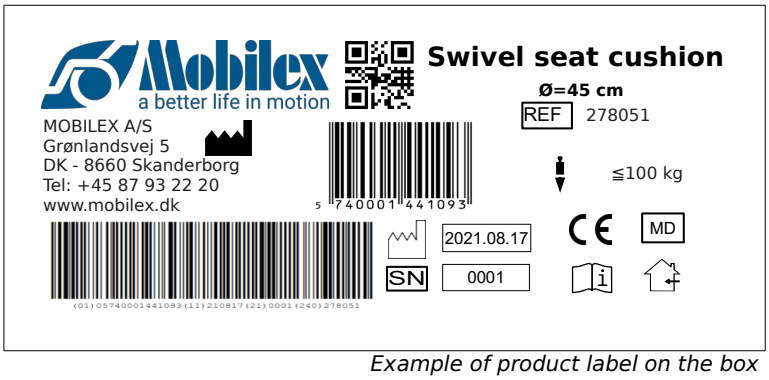Example of product label on the box

5. Care and maintenance
Regular dirt can be removed with standard cleaning products and a sponge or soft rag. Check the specific product information and use only commercial cleaning products which are suitable for cleaning. For disinfection, use only disinfectants certified for use with medical devices available locally. Do not use solvents or abrasives.

Despite the solid construction and the use of resistant materials the product is subject to wear. It is therefore recommended to check the product at regular intervals by a professional service. We recommend the following:

- ✔ Cleaning and disinfection of the product
- ✔ Check the fabric for wear and replace the product if necessary
- ✔ When soiled or stained, the product must be replaced or washed according to the product label
- ✔ The product shall regularly be checked by authorized personnel and always after washing

Recommended service life: Normally 3 years, depending on use.

Storage: We recommend storage in a dry place in temperatures above 0 °C.

Transportation
The device should be transported in a stabilized and secured way with no risk of sudden or dangerous movements during transport.

Reuse
The product is suitable for reuse. Before a new user takes over the product, the device must be checked on a technical level and disinfected according to point 5.

6. Most common problems and solutions

Problems	*Possible causes*	*Solutions*
The swivel seat cushion doesn't turn	The swivel seat cushion has curled up	Uncurl the swivel seat cushion so that it sits flat on the ground

7. Disposal of the product
The product can not be disposed of with household waste and must be brought to the local recycling center.

8. Guarantee terms
Mobilex A/S offers 2 years of warranty for damages caused by errors in production or material errors. The warranty is not valid by not-intended repair or usage of the product. Parts subject to normal wear are not covered by the warranty, unless wear is caused by a manufacturing fault. If product changes are made without our written acceptance the CE-marking and product warranty will not be valid. In the event of product damages covered by this warranty, please notify your dealer or Mobilex A/S directly. The warranty does not cover transport costs and does not include compensation for personal injuries or non-use of the product during repair. The warranty does not cover damages due to the user manual not being followed.

9. Declaration of conformity
The device has passed the conformity assessment in accordance with EU Regulation: MDR 2017/745 regarding medical devices and is CE marked. A declaration of conformity has been issued for the device and is available at mobilex.dk/medias/2021-07/ce_drejesaede__278051_52_278055_56__8274.pdf or via the QR code on the right.

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